

www.drivemedical.com

drive
MEDICAL DESIGN & MANUFACTURING

pneumatic revolving stool



EU Authorized Representative
Drive International, LLC
Globe Lane - Stockport
DUKINFIELD SK16 4RE, GB

www.drivemedical.com


drive
MEDICAL DESIGN & MANUFACTURING

adjustment

To adjust, pull up on handle and move up or down to desired height. Release handle to secure in place.

assembly

1. Remove contents from box.
2. Place base flat on the ground. Insert pneumatic cylinder into base. Make sure the silver end is at the top.
3. Assemble cushion on top of the pneumatic cylinder. Make sure it is secure.
4. Insert castors into the base. Make sure the castors are tightened in the castor holes.

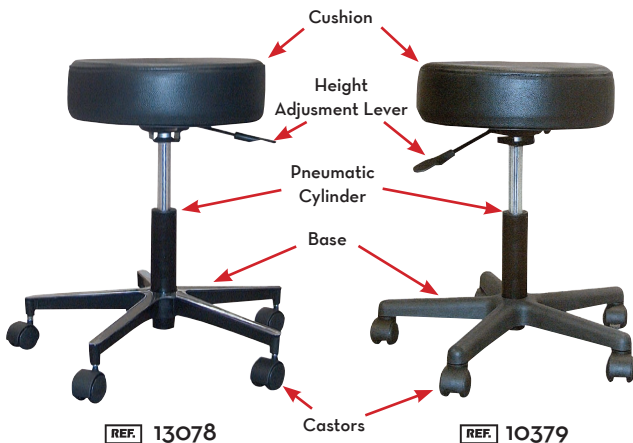

lifetime limited warranty

Your Drive brand product is warranted to be free of defects in materials and workmanship for the lifetime of the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. This Lifetime Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

In the event of a defect covered by this warranty, we will, at our option, repair or replace the device.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.

drive
MEDICAL DESIGN & MANUFACTURING