As metabolism varies among individuals, so does their response to massage. You will soon adjust the time to suit your particular needs and personal comfort. If you have any doubts regarding treatment, we recommend you consult your physician.

## RECOMMENDED MASSAGING TIMES

Generally, massage should last only 15 minutes on any one particular part of the body. Treatment may be applied once or twice daily as desired.

## MAINTENANCE INSTRUCTIONS

- 1) Store massager in a secure place, away from excessive heat or sharp objects. Do not handle the massager by any of the electrical cords and do not wrap the cord around the unit.
- 2) To clean the unit, make sure the massager is unplugged. Allow the unit to cool and then wipe down with a cloth. Do not use any strong chemical cleaners such as benzine.

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- Keep away from children.

- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.

- **Consult your physician should you have any questions regarding therapeutic use of this massager.**
- If you are undergoing treatment for cancer, check with your physician before using.

- This massager should not be used in the following cases: Tuberculosis, benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins, undetermined calf pain, or joint pain.

## CAUTIONS

- Choking hazard. Keep out of reach of children. This is not a toy. Not to be used as a baby or infant pillow. Read instruction booklet before using.

## INSTRUCTIONS FOR USE

1. Place massager around neck.
2. Adjust arms, using one of the 4 positions available, to find a comfortable fit around neck and shoulders.
3. Select desired massage with the easy-to-use control panel. Slide the HEAT button to the ON position to activate heat (function is available with AC adapter only). Turn the VIBRATION dial to choose massage intensity, up to a comfortable level.

**NOTE:** The heat function will work with the AC adapter only, not with the batteries.

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# IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

## DANGER –


To reduce the risk of electric shock:

- 1) **Always unplug this appliance from the electrical outlet immediately after using and before cleaning.**
- 2) **Do not reach for an appliance that has fallen into water.** Unplug immediately.
- 3) Do not use while bathing or in a shower.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

## WARNING –

To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) **An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.**
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) This appliance should not be used by, on, or near children or individuals with certain disabilities.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by Conair.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.

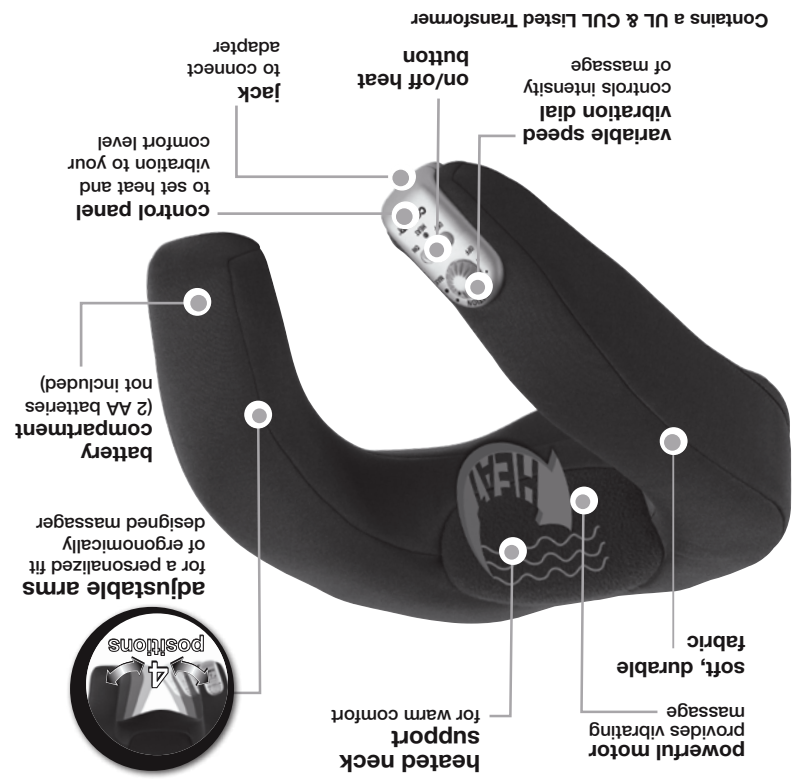**LIMITED ONE YEAR WARRANTY**

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed on back, together with your purchase receipt and \$5.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-366-0937 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO, THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

For information on any Conair® product, call: 1-800-3-CONAIR or visit us on the web at: [www.conair.com](http://www.conair.com)

©2009 Conair Corporation  
 Service Center  
 7475 N. Glen Harbor Blvd.  
 Glendale, AZ 85307  
 09NP14211

NM8X  
 IB-9327

©2009 Conair Corporation  
**CENTRO DE SERVICIO**  
 Conair Corporation  
 Service Department  
 7475 N. Glen Harbor Blvd.  
 Glendale, AZ 85307  
 09NP14211

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